

Store at
-20°C

#33770

SimpleChIP® Human HUWE1 Intron 26 Primers

Small 500 µl (250 PCR reactions)



Support: +1-978-867-2388 (U.S.)
www.cellsignal.com/support

Orders: 877-616-2355 (U.S.)
orders@cellsignal.com


Entrez-Gene ID #10075
UniProt ID #Q7Z6Z7

New 11/18


For Research Use Only. Not For Use In Diagnostic Procedures.

Applications	Species Cross-Reactivity	Primer Anneal/Extension	PCR Product Length
ChIP	H	65°C	94 bp

Description: SimpleChIP® Human HUWE1 Intron 26 Primers contain a mix of forward and reverse PCR primers that are specific to intron 26 of the human HECT, UBA and WWE domain containing 1, E3 ubiquitin protein ligase gene. These primers can be used to amplify DNA that has been isolated using chromatin immunoprecipitations (ChIP). Primers have been optimized for use with SimpleChIP® Universal qPCR Master Mix #88989 and have been tested in conjunction with SimpleChIP® Plus Enzymatic Chromatin IP Kits #9004 and #9005 and ChIP-validated antibodies from Cell Signaling Technology®.

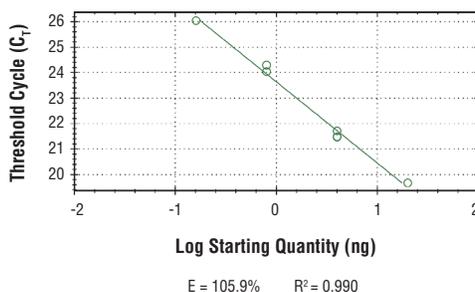

SimpleChIP® Human HUWE1 Intron 26 Primers were tested on DNA isolated from cross-linked cells using the SimpleChIP® Plus Enzymatic Chromatin IP Kit (Magnetic Beads) #9005. Real-time PCR was performed in duplicate on a serial dilution of 2% total input DNA (20 ng, 4 ng, 0.8 ng, and 0.16 ng) using a real-time PCR detection system and SimpleChIP® Universal qPCR Master Mix #88989. The PCR amplification efficiency (E) and correlation coefficient (R²) were calculated based on the corresponding threshold cycle (C_t) of each dilution sample during 40 cycles of real-time PCR (95°C denaturation for 15 sec, 65°C anneal/extension for 60 sec).

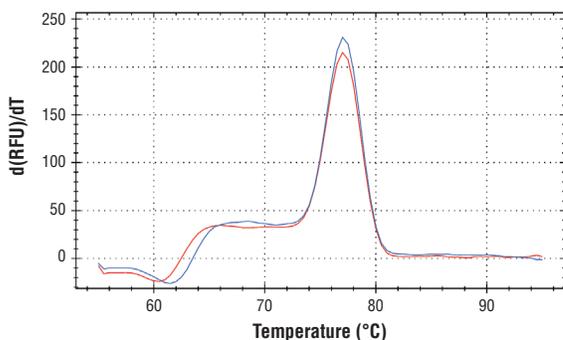

PCR product melting curves were obtained for real-time PCR reactions performed using SimpleChIP® Human HUWE1 Intron 26 Primers. Data is shown for both duplicate PCR reactions using 20 ng of total DNA. The melt curve consists of 80 melt cycles, starting at 55°C with increments of 0.5°C per cycle. Each peak is formed from the degradation of a single PCR product.

Storage: Supplied in nuclease-free water at a concentration of 5 µM (each primer is at a final concentration of 5 µM). Store at -20°C.

Directions for Use:

- Label the appropriate number of PCR tubes or PCR plates compatible with the model of real-time PCR machine to be used. PCR reactions should be performed in duplicate and should include a tube with no DNA to control for contamination, and a serial dilution of a 2% total input chromatin DNA (undiluted, 1:5, 1:25, 1:125), which is used to create a standard curve and determine amplification efficiency.
- Add 2 µl of the appropriate ChIP DNA sample to each tube or well of the PCR plate.
- Prepare a master PCR reaction mix as described below. Add enough reagents for two extra reactions to account for loss of volume. Add 18 µl of the master PCR reaction mix to each PCR reaction tube or well of the PCR plate.

Reagent Volume for 1 PCR Reaction (20 µl)

Nuclease-free H ₂ O	6 µl
5 µM SimpleChIP® Primers	2 µl
SimpleChIP® Universal qPCR Master Mix #88989	10 µl

- Start the following PCR reaction program:
 - Initial Denaturation: 95°C for 3 min.
 - Denaturation: 95°C for 15 sec. c. Anneal and Extension: Primer-specific temp. for 60 sec.
 - Repeat steps b and c for a total of 40 cycles.
- Analyze quantitative PCR results using software provided with the real-time PCR machine.

Background References:

- Orlando, V. (2000) *Trends Biochem Sci* 25, 99-104.
- Kuo, M.H. and Allis, C.D. (1999) *Methods* 19, 425-33.
- Agaloti, T. et al. (2000) *Cell* 103, 667-78.
- Soutoglou, E. and Talianidis, I. (2002) *Science* 295, 1901-4.
- Mikkelsen, T.S. et al. (2007) *Nature* 448, 553-60.
- Lee, T.I. et al. (2006) *Cell* 125, 301-13.
- Weinmann, A.S. and Farnham, P.J. (2002) *Methods* 26, 37-47.
- Wells, J. and Farnham, P.J. (2002) *Methods* 26, 48-56.

Thank you for your recent purchase. If you would like to provide a review visit cellsignal.com/comments.

www.cellsignal.com

© 2018 Cell Signaling Technology, Inc.

SimpleChIP and Cell Signaling Technology are trademarks of Cell Signaling Technology, Inc.

Applications: W—Western IP—Immunoprecipitation IHC—Immunohistochemistry ChIP—Chromatin Immunoprecipitation IF—Immunofluorescence F—Flow cytometry E-P—ELISA-Peptide **Species Cross-Reactivity:** H—human M—mouse R—rat Hm—hamster Mk—monkey Mi—mink C—chicken Dm—D. melanogaster X—Xenopus Z—zebrafish B—bovine Dg—dog Pg—pig S—S. cerevisiae Ce—C. elegans Hr—Horse All—all species expected Species enclosed in parentheses are predicted to react based on 100% homology.