

PathScan® Phospho-p38 MAPK (Thr180/Tyr182) Sandwich ELISA Kit

1 Kit
(96 assays)



Orders ■ 877-616-CELL (2355)
orders@cellsignal.com
Support ■ 877-678-TECH (8324)
info@cellsignal.com
Web ■ www.cellsignal.com

rev. 06/04/19

For Research Use Only. Not For Use In Diagnostic Procedures.

Species Cross-Reactivity: H, M

Description: The PathScan® Phospho-p38 MAPK (Thr180/Tyr182) Sandwich ELISA Kit is a solid phase sandwich enzyme-linked immunosorbent assay (ELISA) that detects endogenous levels of phospho-p38 MAP kinase phosphorylated at Thr180/Tyr182. A phospho-p38 MAP kinase (Thr180/Tyr182) mouse antibody has been coated onto the microwells. After incubation with cell lysates, phospho-p38 MAP kinase (Thr180/Tyr182) protein is captured by the coated antibody. Following extensive washing, a p38 MAP kinase rabbit detection antibody is added to detect the captured phospho-p38 MAP kinase (Thr180/Tyr182). HRP-linked, anti-rabbit antibody is then used to recognize the bound detection antibody. The HRP substrate TMB is added to develop color. The magnitude of the absorbance for this developed color is proportional to the quantity of phospho-p38 MAP kinase (Thr180/Tyr182).

Antibodies in kit are custom formulations specific to kit.

Please visit www.cellsignal.com for a complete listing of recommended companion products.

Specificity/Sensitivity: CST's PathScan® Phospho-p38 MAPK (Thr180/Tyr182) Sandwich ELISA Kit detects endogenous levels of p38 MAP kinase phosphorylated at Thr180/Tyr182 in human and mouse cells. As shown in Figure 1, a significant induction of p38 MAP kinase phosphorylation at Thr180/Tyr182 can be detected in HeLa cells using the PathScan® Phospho-p38 MAPK (Thr180/Tyr182) Sandwich ELISA Kit following treatment with anisomycin or UV irradiation. The level of total p38 MAP kinase remains unchanged as shown by western analysis (Figure 1). This kit detects proteins from the indicated species, as determined through in-house testing, but may also detect homologous proteins from other species.

Background: p38 MAP kinase (MAPK), also called RK (1) or CSBP (2), is the mammalian orthologue of the yeast HOG kinase that participates in a signaling cascade controlling cellular responses to cytokines and stress (1-4). Four isoforms of p38 MAP kinase, p38 α , β , γ (also known as Erk6 or SAPK3), and δ (also known as SAPK4) have been identified. Similar to the SAPK/JNK pathway, p38 MAP kinase is activated by a variety of cellular stresses including osmotic shock, inflammatory cytokines, lipopolysaccharide (LPS), UV light, and growth factors (1-5). MKK3, MKK6, and SEK activate p38 MAP kinase by phosphorylation at Thr180 and Tyr182. Activated p38 MAP kinase has been shown to phosphorylate and activate MAPKAP kinase 2 (3) and to phosphorylate the transcription factors ATF-2 (5), Max (6), and MEF2 (5-8). SB203580 (4-(4-fluorophenyl)-2-(4-methylsulfinylphenyl)-5-(4-pyridyl)-imidazole) is a selective inhibitor of p38 MAPK. This compound inhibits the activation of MAPKAP-2 by p38 MAPK and subsequent phosphorylation of HSP27 (9).

Product Includes	Item #	Kit Quantity	Color	Storage Temp
Phospho-p38 MAPK (T180/Y182) MmAb Coated Microwells*	81256	96 tests		4°C
p38 MAPK Detection Rabbit mAb	32613	1 each	Green (Lyophilized)	4°C
Anti-rabbit IgG, HRP-linked Antibody (ELISA Formulated)	13272	1 each	Red (Lyophilized)	4°C
Detection Ab Diluent 3	14632	11 ml	Green	4°C
HRP Diluent	13515	11 ml	Red	4°C
TMB Substrate	7004	11 ml		4°C
STOP Solution	7002	11 ml		4°C
Sealing Tape	54503	2 sheets		4°C
ELISA Wash Buffer (20X)	9801	25 ml		4°C
ELISA Sample Diluent	11083	25 ml	Blue	4°C
Cell Lysis Buffer (10X)	9803	15 ml		-20°C

*12 8-well modules – Each module is designed to break apart for 8 tests.

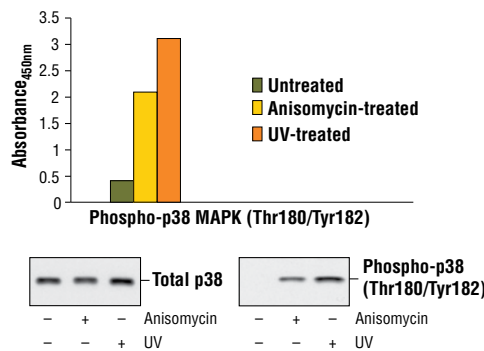


Figure 1. Treatment of HeLa cells with anisomycin or UV irradiation stimulates phosphorylation of p38 MAP kinase at Thr180/Tyr182 as detected by the PathScan® Phospho-p38 MAPK (Thr180/Tyr182) Sandwich ELISA Kit #7946, but does not affect the level of total p38 MAP kinase. HeLa cells were treated with 5 μ g/ml anisomycin for 30 minutes at 37°C or 45 mJ/cm² UV irradiation followed by a 30 minute recovery period at 37°C. The absorbance readings at 450 nm are shown in the top figure, while the western blots using p38 MAPK Antibody #9212 (left) or Phospho-p38 MAPK (Thr180/Tyr182) Antibody #9211 (right) are shown in the bottom figure.

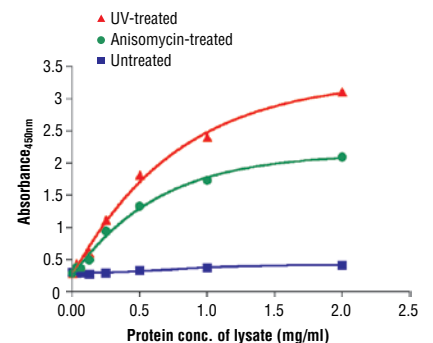


Figure 2. The relationship between lysate protein concentration from untreated, anisomycin-treated or UV-irradiated HeLa cells and the absorbance at 450 nm is shown.

SB203580 inhibits p38 MAPK catalytic activity by binding to the ATP binding pocket, but does not inhibit phosphorylation of p38 MAPK by upstream kinases (10).

PathScan® Sandwich ELISA Protocol (for kits with Lyophilized Antibodies)

A Solutions and Reagents

NOTE: Prepare solutions with purified water.

- Microwell strips:** Bring all to room temperature before use.
- Detection Antibody:** Supplied lyophilized as a green colored cake or powder. Add 1.0 ml of Detection Antibody Diluent (green solution) to yield a concentrated stock solution. Incubate at room temperature for 5 min with occasional gentle mixing to fully reconstitute. To make the final working solution, add the full 1.0 ml volume of reconstituted Detection Antibody to 10.0 ml of Detection Antibody Diluent in a clean tube and gently mix. Unused working solution may be stored for 4 weeks at 4°C.
- HRP-Linked Antibody*:** Supplied lyophilized as a red colored cake or powder. Add 1.0 ml of HRP Diluent (red solution) to yield a concentrated stock solution. Incubate at room temperature for 5 min with occasional gentle mixing to fully reconstitute. To make the final working solution, add the full 1.0 ml volume of reconstituted HRP-Linked Antibody to 10.0 ml of HRP Diluent in a clean tube and gently mix. Unused working solution may be stored for 4 weeks at 4°C.
- Detection Antibody Diluent:** Green colored diluent for reconstitution and dilution of the detection antibody (11 ml provided).
- HRP Diluent:** Red colored diluent for reconstitution and dilution of the HRP-Linked Antibody (11 ml provided).
- Sample Diluent:** Blue colored diluent provided for dilution of cell lysates.
- 1X Wash Buffer:** Prepare by diluting 20X Wash Buffer (included in each PathScan® Sandwich ELISA Kit) in purified water.
- Cell Lysis Buffer:** 10X Cell Lysis Buffer #9803 or 1X Cell Lysis Buffer #7018: This buffer can be stored at 4°C for short-term use (1–2 weeks). Recommended: Add 1 mM phenylmethylsulfonyl fluoride (PMSF) immediately before use.
- TMB Substrate** (#7004).
- STOP Solution** (#7002).

*Note: Some PathScan® ELISA Kits may include HRP-Linked Streptavidin in place of HRP-Linked Antibody.

B Preparing Cell Lysates

For adherent cells.

- Aspirate media when the culture reaches 80–90% confluence. Treat cells by adding fresh media containing regulator for desired time.
- Remove media and rinse cells once with ice-cold 1X PBS.
- Remove PBS and add 0.5 ml ice-cold 1X Cell Lysis Buffer plus 1 mM PMSF to each plate (10 cm diameter) and incubate the plate on ice for 5 min.
- Scrape cells off the plate and transfer to an appropriate tube. Keep on ice.
- Sonicate lysates on ice.
- Microcentrifuge for 10 min (14,000 rpm) at 4°C and transfer the supernatant to a new tube. The supernatant is the cell lysate. Store at –80°C in single-use aliquots.

For suspension cells

- Remove media by low speed centrifugation (~1200 rpm) when the culture reaches 0.5–1.0 × 10⁶ viable cells/ml. Treat cells by adding fresh media containing regulator for desired time.
- Collect cells by low speed centrifugation (~1200 rpm) and wash once with 5–10 ml ice-cold 1X PBS.
- Cells harvested from 50 ml of growth media can be lysed in 2.0 ml of 1X Cell Lysis Buffer plus 1 mM PMSF.
- Sonicate lysates on ice.
- Microcentrifuge for 10 min (14,000 rpm) at 4°C and transfer the supernatant to a new tube. The supernatant is the cell lysate. Store at –80°C in single-use aliquots.

C Test Procedure

- After the microwell strips have reached room temperature, break off the required number of microwells. Place the microwells in the strip holder. Unused microwells must be resealed and stored at 4°C immediately.
- Cell lysates can be undiluted or diluted with Sample Diluent (supplied in each PathScan® Sandwich ELISA Kit, blue color). Individual datasheets for each kit provide a sensitivity curve that serves as a reference for selection of an appropriate starting lysate concentration. The sensitivity curve shows typical kit assay results across a range of lysate concentration points.
- Add 100 µl of each undiluted or diluted cell lysate to the appropriate well. Seal with tape and press firmly onto top of microwells. Incubate the plate for 2 hr at 37°C. Alternatively, the plate can be incubated overnight at 4°C.
- Gently remove the tape and wash wells:
 - Discard plate contents into a receptacle.
 - Wash 4 times with 1X Wash Buffer, 200 µl each time for each well.
 - For each wash, strike plates on fresh towels hard enough to remove the residual solution in each well, but do not allow wells to completely dry at any time.
 - Clean the underside of all wells with a lint-free tissue.
- Add 100 µl of reconstituted Detection Antibody (green color) to each well (refer to Section A, Step 2). Seal with tape and incubate the plate at 37°C for 1 hr.
- Repeat wash procedure (Section C, Step 4).
- Add 100 µl of reconstituted HRP-Linked secondary antibody (red color) to each well (refer to Section A, Step 3). Seal with tape and incubate the plate for 30 min at 37°C.
- Repeat wash procedure (Section C, Step 4).
- Add 100 µl of TMB Substrate to each well. Seal with tape and incubate the plate for 10 min at 37°C or 30 min at 25°C.
- Add 100 µl of STOP Solution to each well. Shake gently for a few seconds.

NOTE: Initial color of positive reaction is blue, which changes to yellow upon addition of STOP Solution.

- Read results.
 - Visual Determination:** Read within 30 min after adding STOP Solution.
 - Spectrophotometric Determination:** Wipe underside of wells with a lint-free tissue. Read absorbance at 450 nm within 30 min after adding STOP Solution.

Background References:

- Rouse, J. et al. (1994) *Cell* 78, 1027–1037.
- Han, J. et al. (1994) *Science* 265, 808–811.
- Lee, J.C. et al. (1994) *Nature* 372, 739–746.
- Freshney, N.W. et al. (1994) *Cell* 78, 1039–1049.
- Raingaud, J. et al. (1995) *J. Biol. Chem.* 270, 7420–7426.
- Zervos, A.S. et al. (1995) *Proc. Natl. Acad. Sci. USA* 92, 10531–10534.
- Zhao, M. et al. (1999) *Mol. Cell. Biol.* 19, 21–30.
- Yang, S.H. et al. (1999) *Mol. Cell. Biol.* 19, 4028–4038.
- Cuenda, A. et al. (1995) *FEBS Lett* 364, 229–33.
- Kumar, S. et al. (1999) *Biochem Biophys Res Commun* 263, 825–31.